
MANUAL DE INSTRUCCIONES

Mx Onda

IRRIGADOR DENTAL

MODELO: **MX-HD2424**

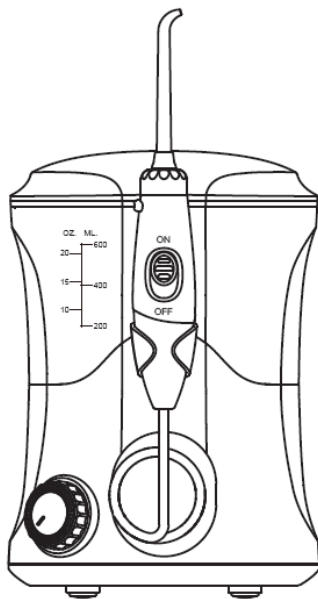CE

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
MODEL: **MX-HD2424**

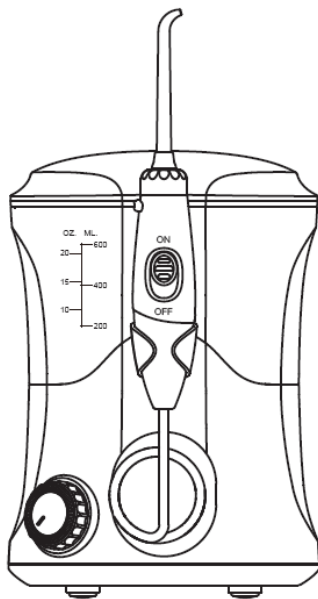CE

BEFORE USE READ THIS MANUAL CAREFULLY

ATTENTION

THIS APPLIANCE WORKS WITH A VOLTAGE OF 230 V, TO AVOID POSSIBLE ELECTRIC SHOCK DO NOT TRY TO OPEN IT. This equipment does not have any part that can be repaired by yourself, in case of a malfunction do not handle it, go to an Official Technical Service MX ONDA.

Once the useful life of this electrical product has been exhausted, do not throw it into the household waste. Deposit it in the appropriate container of a clean point in your locality, for later recycling.

This product complies with the European RoHS Directive (2011/65 / EU), on the restriction of the use of certain dangerous substances in electrical and electronic equipment.

INTRODUCTION

MX ONDA thanks you for the deference you have had when purchasing this product. The use of this dental irrigator will help keep your teeth and gums in perfect health.

Regarding dental care, many people tend to measure the health of their teeth by observing only the whiteness of the teeth, ignoring the importance of the health of the gums.

The failure to clean only the surface of the teeth, and not the gums or the interdental space, apart from causing halitosis, causes the accumulation of plaque and tartar on the teeth. A continued poor dental hygiene can cause periodontal diseases such as gingivitis or periodontitis.

Taking into account the importance of proper dental cleaning, the continued use of this device after each meal, is a first step to achieve a perfect dental hygiene.

With this dental irrigator along with the supplied nozzles, you can effectively clean the teeth, the interdental space, the gum line, the tongue and, if necessary, the orthodontics.

By orienting the nozzle you can adjust the direction of water output, thus accessing areas that are not easily reached with the traditional brush. Thanks to the pressure with which the water is projected, while cleaning the teeth will produce a pleasant massage on the gums, helping to improve blood circulation and preventing dental diseases.

BEFORE CONNECTING THE DENTAL IRRIGATOR

When using electrical appliances, basic safety precautions should always be followed to reduce the risk of fire, electric shock and / or bodily injury. These measures include the following:

- Read these instructions carefully and keep them for later use.
- Check that the voltage of your outlet is 230V ~ 50Hz.
- To protect against the risk of electric shock and personal injury, do not insert the cord, plug or base of the appliance in water or any other liquid.
- To reduce the risk of fire or electric shock, do not remove the cover of the appliance.
- After removing the packaging, check that the appliance is not damaged. If in doubt, use it; Contact the technical assistance staff.
- Packaging materials (bags, etc.) should not be left within the reach of children, as they are a potential source of danger.
- If you are following any kind of oral treatment or suffer from periodontal disease, consult your usual dentist before using this device.

SECURITY MEASURES

- This product is not a toy, children should not play with this device. The cleaning and maintenance to be performed by the user can not be performed by children, unless they are over 8 years old and are supervised.


- This device can not be used by children (under 8 years old) or people with reduced physical, sensory or mental abilities or without experience and knowledge, unless they have received appropriate instructions on the use of the device and are supervised by an adult or responsible person.
- When handling this appliance, never do so with wet hands or feet or inside the bathtub, shower, etc.
- Do not insert the device in water or any other liquid.
- This device is only designed for domestic use. MX ONDA will not accept any responsibility and the guarantee will be invalid in case of improper use or non-compliance with the instructions.
- Do not leave the device exposed to the elements (rain, sun, etc.).
- Do not use the appliance outdoors or in places where aerosol products (sprays) or oxygen are used.
- Do not use any accessory that is not recommended by the manufacturer.
- Do not use the appliance if the nozzles are damaged or plugged.
- Do not obstruct the air inlets of the irrigator.
- While the irrigator is connected, do not neglect your vigilance and keep it out of reach of children.
- If the dental irrigator is used by several people, for reasons of hygiene it is recommended that each person use their own mouthpiece.
- To unplug the cable from the network, gently pull the plug from it. Never do it by pulling the cable.
- Do not twist the tube or wrap it around the device.
- Do not apply the water jet directly under the tongue, inside the ears or in the nostrils. High pressure can damage sensitive areas.
- This product is only for dental buco hygiene. Please do not apply it to other areas.
- Unplug the irrigator whenever it is not in use or absent from the place (even for a brief moment), when it is connected to the mains and there is a power failure, or "jump" a fuse or automatic and before perform maintenance cleaning itself.

- Use the appropriate nozzle for the area and type of cleaning that you want to carry out, choosing an inadequate nozzle can cause damage to the applied areas. Do not use the irrigator to clean the ears or other sensitive areas.
- In case of failure, anomalies, if the device has been dropped, the plug or other parts are damaged or the network cable is damaged, do not try to repair it on your own. Contact the nearest MX ONDA Official Technical Service.


TIPS

The irrigator includes 7 tips; 3 standard tips, 1 periodontal tip, 1 orthodontic mouthpiece tip, 1 mouthpiece for cleaning dental plaque and 1 tip for cleaning the tongue.

The tips have been designed to obtain the best results in cleaning the area for which they were made as well as their easy and comfortable use.

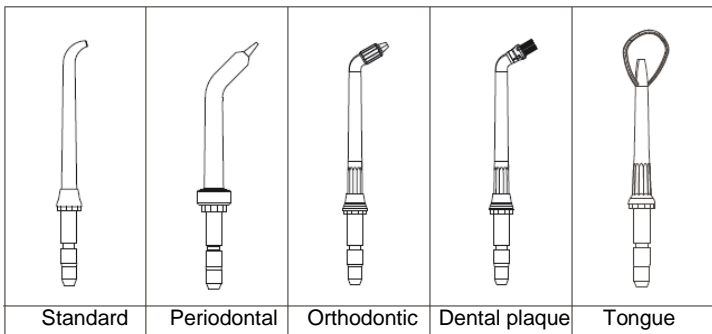CONTROLS AND ELEMENTS

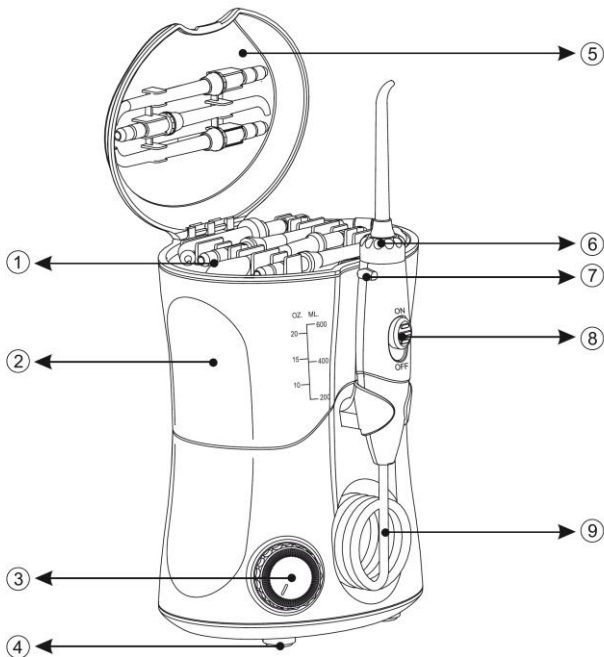
1. Reservoir cover and tips compartment
2. Reservoir
3. On / off switch and pressure regulation
4. Anti-slip feet
5. Tips compartment cover
6. 360 ° tip rotation
7. Release button of tip
8. Blocking the water outlet
9. Self-coiling tube

USE OF THE ORAL IRRIGATOR

Filling the reservoir

1. Check that the irrigator is unplugged from the mains socket.
2. Remove the reservoir **2** from the base and remove the lid **1** from the reservoir.
3. Fill the reservoir without exceeding the maximum capacity of the same, you can use water from time, hot water (without the temperature of the same exceeds 40 °C) or mixed with some type of mouthwash or liquid mouthwash.
4. Replace the cover in the reservoir, and place the reservoir in the base, pressing it lightly until it is correctly fitted.

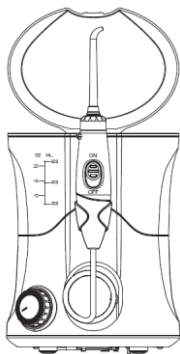Inserting and removing the tips

To insert the tip in the handle, match the axis of the handle with the hole in the handle, then push the tip towards the handle (figure **A**), you will notice that the rotation ring of the tip **6** slides down, and You will hear a "click" signal that the tip is properly retained.

To remove the tip, press the release button **7**, then pull the tip upwards (figure **B**).

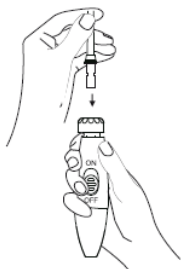
Figure A (Inserting)


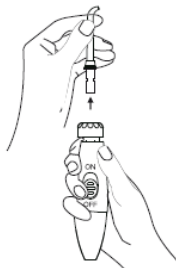
Figure B (Removing)

Cover with tips compartment

The cover **1** of the reservoir includes a compartment for storing the tips. Inside the compartment is engraved in relief the position of each tip, to facilitate its placement.


360 ° Tip rotation

To facilitate cleaning and access to any area of the mouth, this irrigator incorporates a device that allows you to rotate the tip 360°. Turn the rotation ring of the tip **6** and orient it according to the type or cleaning zone that you wish to make.

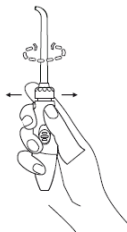Blocking the water outlet

When you need to stop the water momentarily, to change the position of the tip or to rinse, slide the button for blocking the water outlet **8** located on the handle to the **OFF** position, to allow the exit of the water slide this button towards **ON** position.

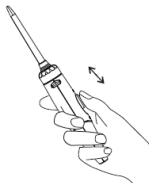
On / off and pressure regulation

To turn **ON** the irrigator, turn control **3** to the right (you will hear a "click").

To increase the water pressure turn control **3** to the right (clockwise) or turn it to the left (counter clockwise) to decrease the water pressure.

To turn **OFF** the dental irrigator, turn control **3** to the left until you hear a "click".

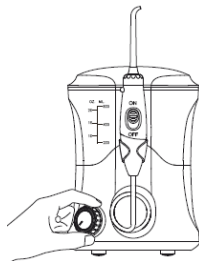You can adjust the pressure depending on the type of cleaning you want to perform and also the level of tooth and gum sensitivity, which varies or is different in each person.

Note:


We recommend that you start with a low water pressure level, then increase the pressure level according to your preference.

Function of each tip

The irrigator includes 5 different tips; standard mouthpiece, periodontal mouthpiece, for orthodontics, for cleaning the tongue and for cleaning dental plaque. The tips have been designed to obtain the best results in the cleaning of the area for which it was made as well as its easy and comfortable use.

Standard tip


Use this mouthpiece to clean the interdental cavities (space between the teeth), you can also clean the line between the gum and the teeth.

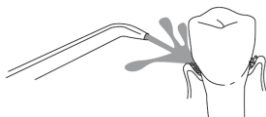*Periodontal tip*

The periodontal mouthpiece has a rubber tip, so that it can be inserted between the gum and the teeth and does not cause pain or discomfort, as the rubber adapts to the small space. The rubber tip must be at 45 ° to the tooth.

*Orthodontic tip*

This nozzle has a small brush in the tip that together with the water jet, will perform a deep cleaning of the elements that make up the orthodontics, eliminating the remains deposited between the teeth and the orthodontic element.

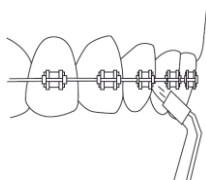
Tongue tip

How to clean the tongue does not need a direct stream of water, this nozzle is shaped like a spoon in order to diffuse the jet of water in a fan, thus not only increases the surface area to clean, it adapts the water jet to clean the tongue.

Tip for cleaning dental plaque

Direct the jet at a 90 ° angle to the gum line, close the lips slightly. Begin at the posterior molar area and towards the teeth, place the tip near the teeth so that the bristles gently touch the teeth.


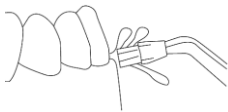
Note:


Until you become familiar with the irrigator, for any type of cleaning; dental, periodontal, orthodontics or for the tongue, we recommend that you start with a low level of water pressure, then increase the pressure level according to your preference.

USE OF ORAL IRRIGATOR

1. Check that the irrigator is off and the control **3** is in the OFF position.
2. Once the reservoir is filled with water or solution and with the appropriate tip to the type of cleaning to be done attached to the handle, connect the irrigator to a mains power socket.
3. Tilt your head over the sink and keep it at an approximate 45 degree angle. Insert the outlet of the tip into the mouth.
4. Slide the water outlet lock button **8** to the **OFF** position and turn on the irrigator by turning control **3** clockwise until a "click" is heard.

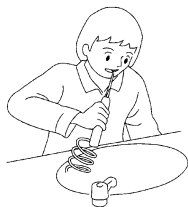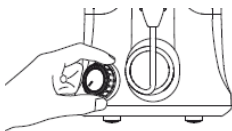5. Slide the lock button on the water outlet **8** to the **ON** position and the water or the solution introduced in the reservoir **2**, will begin to exit through the tip.
6. During cleaning, you can increase or decrease the water flow pressure by means of control **3**, or turn the tip by means of the rotation ring of tip **6** to direct the water jet towards the desired zone.
7. When you need to stop the water coming out, to remove the tip from the mouth to rinse, rest, etc. Slide the locking button of the water outlet **8** to the **OFF** position. Never remove the tip from the mouth if you have not turned off the irrigator before or stopped the water flow by pressing button **8**.
8. When cleaning is complete, slide the water outlet lock button **8** to the **OFF** position, remove the mouthpiece from the mouth, and then turn the appliance off by turning control **3** clockwise until a "click" is heard.
9. Unplug the irrigator from the mains socket, remove reservoir **2** and empty it in the sink.
10. Replace the reservoir **2** in the base, point the tip towards the sink and turn on the irrigator by turning control **3** to the right, this will empty the pump and the flexible tube. When the water stops flowing, turn the appliance off again by turning control **3** to the left until you hear a "click".
11. When emptying the pump and the flexible tube, it will prevent the formation of bacteria inside the pump and the flexible tube due to the accumulated water.
12. Place the handle in its holder, taking care that the self-coiling tube is correctly lodged in its housing and that it is not bent or forced.
13. Remove the tip and store it in the compartment located on the reservoir lid and **9**. Unplug the irrigator from the mains socket.

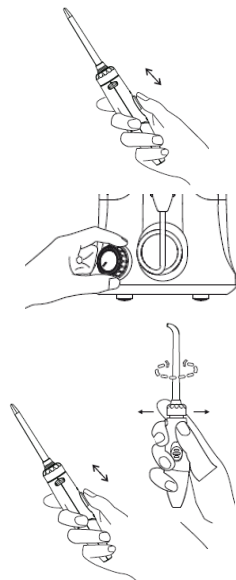Very important:

- Water splashes may occur during the use of this device. Start by first selecting the lowest water flow pressure and always keep the head tilted over the sink. This way you will avoid an excess of water spill out of the sink.
- To avoid overheating the pump, never run the irrigator without water or solution for a long time.
- When possible, do not use this device continuously for more than 6 minutes at a time. Take breaks of at least 1 minute.
- Unplug the appliance from the mains socket before filling or emptying the reservoir.

TIMER

Normally the time required for a full and deep oral cleaning is approximately 3 minutes. The dental irrigator incorporates a 3-minute timer and consists in automatically turning it off after that time.

How does the timer work?

If during the cleaning process the irrigator is not turned off by control 3, after 3 minutes the irrigator will automatically turn off. If you wish to continue using the irrigator, you must turn off the irrigator by turning control **3** to the left until you hear a "click", then turn it on again by turning this control clockwise.

BASIC DENTAL CARE

If a slight bleeding of the gums occurs during the use of this product, it means that you may have a periodontal problem.

Continue using this irrigator for a couple of weeks, the bleeding of the gums should disappear, if not stop using this product and consult your usual dentist.


The pressure that this device provides is approximately 1700 impulses per minute, being able to effectively clean the interdental cavities, massage the gums, eliminate halitosis and prevent possible dental diseases.

If the user has implanted an orthodontic appliance, it is recommended that he perform the cleaning in a way, even more conscientiously, than the rest of the users. In this way it will prevent the traces of particles from depositing between the teeth.

For proper dental care you should use this device, after each meal, you must also visit your regular dentist on a regular basis.

Cleaning the dental surface

- Turn the outlet of the tip in such a way that the jet of water forms an angle of 90 degrees with respect to the teeth, see the following illustration.
- Change the position of tip to clean all the teeth.
Clean the teeth with great care, especially paying attention to the interdental spaces.

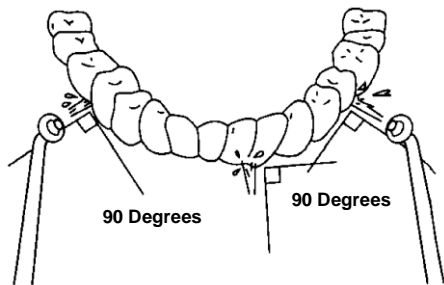Cleaning the inside of the denture

- Orient the tip so that the water outlet is directed towards the inside of the teeth. Make sure that the jet of water forms an angle of 90 degrees with respect to the tooth.
- The cleaning of the inside of the denture, both upper and lower (Fig.1 and Fig.2), is done in a similar way. Always keep your mouth slightly open, to facilitate the exit of water from your oral cavity. As indicated above, the head should always be slightly tilted above the sink.

Fig.1


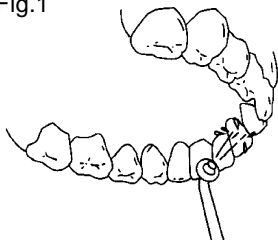
Fig.2

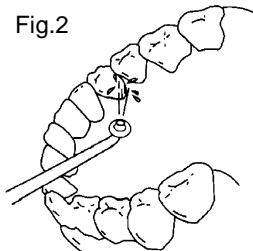CLEANING AND MAINTENANCE

- Unplug the irrigator from the mains socket.
- Clean the surfaces of the irrigator with a damp cloth and dry them with a dry one. Do not use alcohol or abrasive products.
- Do not insert the device in water or any other liquid. Do not wrap the tube around the appliance.
- Store the tips in their compartment or in an area free of dust or dirt.

¡REMEMBER!

Never handle or use the irrigator with wet hands, or while showering or bathing. In this way the possibility of electric shocks is avoided.

TECHNICAL SPECIFICATIONS

Supply voltage	100 - 240 V ~ 50/60 Hz
Consumption	18 W
Capacity	0.6 liters
Length of use of the tube	80 cm approx.
Pressure	30 - 90 PSI
Pulses per minute	1700
Dimensions	140 x 125 x 200 mm (l, a, height)
Weight	750 g

Specifications subject to change without notice.

This product complies with the European Directives 2004/108/EC (EMC – Electromagnetic Compatibility) and the 2006/95/EC (LVD – Low Voltage Directive) relative to Low Tension Security.

“MX, MX ONDA” and its logos are trademarks of MX ONDA, S.A.

Telephone of information and Technical Service: **+34 902 551 501**

MX ONDA, S.A.

Isla de Java, 37

28034 – MADRID (SPAIN)

Manual de instrucciones
Instruction manual
Manual de instruções
Benutzerhandbuch
Manuel utilisateur
Manuale di istruzioni

E-MAIL: mxsat@mxonda.es

<http://www.mxonda.es>